

FULL SIZE PORTABLE SCOOTER



Seat Color  
Black Vinyl

Body Colors  
Blue  
Red

**LiteRider**<sup>GOLDEN</sup>  
3 wheel GL-111  
4 wheel GL-141

Enjoy the best of both worlds with the LiteRider™ scooter! The LiteRider is a full-size scooter with all of the benefits of a compact, portable scooter! The LiteRider comes standard with a Delta Tiller and a large, 17" x 16" stadium style seat and adjustable armrests. With a 300 lb. weight capacity, the LiteRider can be quickly and wirelessly disassembled for transport or storage.



Model	GL111	GL141
Medicare Code	N/A	N/A
Weight Capacity	300 lbs	300 lbs
Drive Wheels	Rear	Rear
Maximum Speed	5 mph	5 mph
Operating Range up to <sup>1</sup>	12 miles	12 miles
Ground Clearance to center deck	3.75"	3.75"
Ground Clearance to transaxle/motor	2.5"	2.5"
Turning Radius	39"	48"
Type Batteries	2-22AH	2-22AH
Colors	Red, Blue	Red, Blue
Freewheel Mode	Yes	Yes
Electronic Speed Control	Yes	Yes
Electro-Mechanical Brakes	Yes	Yes
Charger	Off Board	Off Board
Controller	90 Amp Dynamic R Series	90 Amp Dynamic R Series

Dimensions <sup>3</sup>		
Length	41.5"	43.25"
Width	23.25"	22"
Height (ground to top of back of seat)	37" - 39"	37" - 39"
Ground to Top of Seat	23.5" - 25.5"	23.5" - 25.5"
Ground to Top of Deck	6"	6"
Deck to Top of Seat	18" - 20"	18" - 20"
Number of Seat Height Positions	3	3
Size of Increments	1"	1"
Front Axle to Rear Axle	30"	30"
Weight of Unit: (assembled)		
Front Half	44 lbs	52 lbs
Rear Half	28 lbs	28 lbs
Seat w/Arms	22 lbs	22 lbs
Batteries	29 lbs/pack	29 lbs/pack
Tires:		
Front	Solid/Low Profile	Solid/Low Profile
Rear	9"	9"
Rear Anti-Tippers	2"	2"
Standard Seat <sup>3</sup>		
Back Height - no headrest	Stadium Style	Stadium Style
Width x Depth	16"	16"
Color	17" X 16"	17" X 16"
	Black Vinyl	Black Vinyl

Warranty:		
Frame	Lifetime	Lifetime
Drive Train	2 Years	2 Years
Electronics <sup>2</sup>	2 Years	2 Years
Batteries	1 Year	1 Year

<sup>1</sup>. Battery range tested at 200 lbs., and will vary due to rider weight, drive surface and terrain.

<sup>2</sup>. Electronics warranty excludes batteries.

## Features

- Delta Tiller standard
- 300 lb. weight capacity
- Rugged anti-rust steel frame tubing
- Rugged plastic floor mat for optimum durability
- Disassembles quickly and easily

## Enjoy quick and wireless disassembly!




Pick up the seat



Pull out the battery pack



Remove the baskets from under the seat and from the tiller (not shown)

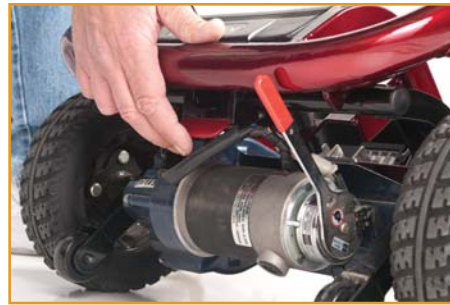

To separate the rear transaxle from the frame, squeeze the release lever located on the underside of the frame




While squeezing the release lever, pull up on the rear frame



Disassembles into six pieces for transport and storage! No wires to connect or disconnect!



Bright adjustable headlight

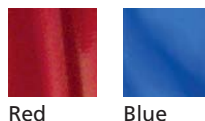

Red

Blue

Literature is current at time of printing. Golden Technologies reserves the right to make changes to the product or literature at any time.